

# ARIADNE'S THREAD®

## Repositionable Breast Localisation Needle



*With its retractable hook and the unique flexibility of its guiding thread, Ariadne's Thread® is an effective tool in pinpointing microcalcification and mammary lesions.*

### ARIADNE'S THREAD® & ARIADNE'S THREAD II®


#### ► A straightforward procedure...

- **Wing grip and thumbbutton design** allows single handed use,
- **Retractable hook** allows repositioning,
- **Indicator on wing** shows direction of hook.



#### ► ...for accurate location...

- **Hook release at the needle tip** for an accurate pinpointing,
- **Break indicators & silicone stopper** for accurate depth location,
- **Echogenic Needle Tip** for ultrasound use,
- **Removable silicone stopper** for partial hook release prior to confirming position.

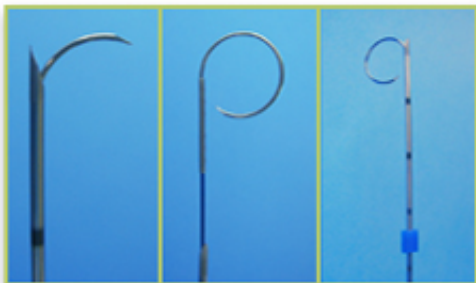

#### ► ...and efficient surgery...

- **8mm Loop Hook** for easier penetration and stability,
- **Metal or flexible Polyamide guide thread** to reduce movement,
- **Silicone limiter** for partial hook release prior to confirming position.

#### ► References ARIADNE'S THREAD® Products

*Ariadne's Thread is available in three models:*

- **Ariadne's Thread®** with polyamid thread
- **Ariadne's Thread®** with metal thread
- **Ariadne's Thread II®** with metal thread

CAT. NUMBER	ARIADNE'S THREAD	LENGTH	THREAD
FA 20/7	Ariadne'sthread L70	70mm	Polyamid
FA 20/10	Ariadne'sthread L100	100mm	Polyamid
FA 20/7	Ariadne'sthread L70	70mm	Metal
FA 20/10	Ariadne'sthread L100	100mm	Metal
FA 20/135	Ariadne'sthread L135	135mm	Metal