

Laurane™ Bone Biopsy Kit

Bone Biopsy Kit	KBD 11/10	KBD 11/15
Aenesthesia needle	L90mm	L150mm
Stylet with Trocar Tip	2.4mm x 140mm	2.4mm x 190mm
Cannula	11G x 100mm	11G x 150mm
Biopsy needle with sharp trephined tip	2.35mm x 193mm	2.35mm x 233mm
Needle obturator	1.7mm x 234mm	1.7mm x 274mm
Coaxial Bone Biopsy Kit	KBC 11/10	KBC 11/15
Aenesthesia needle	L90mm	L150mm
Guiding Stylet	0.9mm x 250mm	0.9mm x 300mm
Coaxial Needle	2.3mm x 138mm	2.3mm x 188mm
Cannula	11G x 100mm	11G x 150mm
Biopsy needle with sharp trephined tip	2.35mm x 193mm	2.35mm x 233mm
Needle obturator	1.7mm x 234mm	1.7mm x 274mm
Hard Bone Biopsy Kit	KBP 11/10	KBP 11/15
Aenesthesia needle	L90mm	L150mm
Stiffening Stylet	2.3mm x 250mm	2.3mm x 300mm
Sharp Tip Cannula	11G x 100mm	11G x 150mm
Drill	2.4mm x 137mm	2.4mm x 187mm
Biopsy needle with sharp trephined tip	2.35mm x 193mm	2.35mm x 233mm
Needle obturator	1.7mm x 234mm	1.7mm x 274mm
Laurane™ Perforating Cannula Drill	KTP 11/10	KTP 11/15
Stiffening Stylet	2.3mm x 250mm	2.3mm x 300mm
Sharp Tip Cannula	11G x 100mm	11G x 150mm
Drill	2.4mm x 137mm	2.4mm x 187mm
Biopsy Needle	CD 11/10	CD 11/15
Biopsy needle with sharp trephined tip	2.35mm x 193mm	2.35mm x 233mm
Needle obturator	1.7mm x 234mm	1.7mm x 274mm
Long Drill	ML 11/10	ML 11/15
Long drill with handle	2.4mm x 186mm	2.4mm x 236mm

Products	Order No.	Diameter	Length
Bone Biopsy Kit	KBD 11/10	11G	100mm
Bone Biopsy Kit	KBD 11/15	11G	150mm
Coaxial Bone Biopsy Kit	KBC 11/10	11G	100mm
Coaxial Bone Biopsy Kit	KBC 11/15	11G	150mm
Hard Bone Biopsy Kit	KBP 11/10	11G	100mm
Hard Bone Biopsy Kit	KBP 11/15	11G	150mm
Laurane™ Perforating Cannula Drill	KTP 11/10	11G	100mm
Laurane™ Perforating Cannula Drill	KTP 11/15	11G	150mm
Biopsy needle	CD 11/10	2.35mm	193mm
Biopsy needle	CD 11/15	2.35mm	233mm
Long drill	ML 11/10	2.4mm	186mm
Long drill	ML 11/15	2.4mm	236mm

Prof. J.D. Laredo*

Dr. B. Hamze*



*Department of Radiology - Lariboisiere Hospital, Paris-France.




Many thanks to Professor J.D. Laredo and Dr. B. Hamze for their most helpful contribution to the development of the Laurane™ Bone Biopsy Kit.

<http://www.lauranemedical.com>



Optional Accessories

Laurane™ Perforating Cannula Drill: KTP

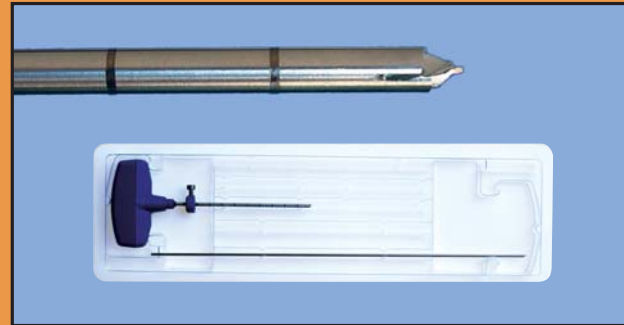

May be used in addition with KBD or KBC when bone appeared to be hard bone.

- ◆ Stiffening Stylet
- ◆ Sharp Tip Cannula
- ◆ Drill

Long Drill: ML



May be used in addition with KBP or KTP for deep lesions in hard bone.

Biopsy Needle: CD




May be used in addition with KBD/KBC/ KBP when more than two core samples are required.

- ◆ Biopsy Needle With Sharp Trepined Tip
- ◆ Needle Obturator

Three Kits to Ensure a Perfect Procedure

Bone Biopsy Kit: KBD


- ◆ Anesthesia Needle
- ◆ Stylet With Trocar Tip
- ◆ Cannula
- ◆ Biopsy Needle With Sharp Trepined Tip
- ◆ Needle Obturator

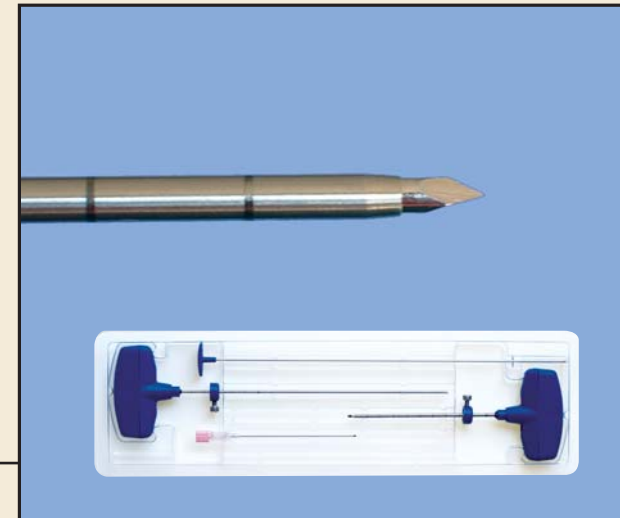

The Bone Biopsy Kit is used to biopsy superficial bones.

Coaxial Bone Biopsy Kit: KBC


- ◆ Anesthesia Needle
- ◆ Guiding Stylet
- ◆ Coaxial Needle
- ◆ Cannula
- ◆ Biopsy Needle With Sharp Trepined Tip
- ◆ Needle Obturator

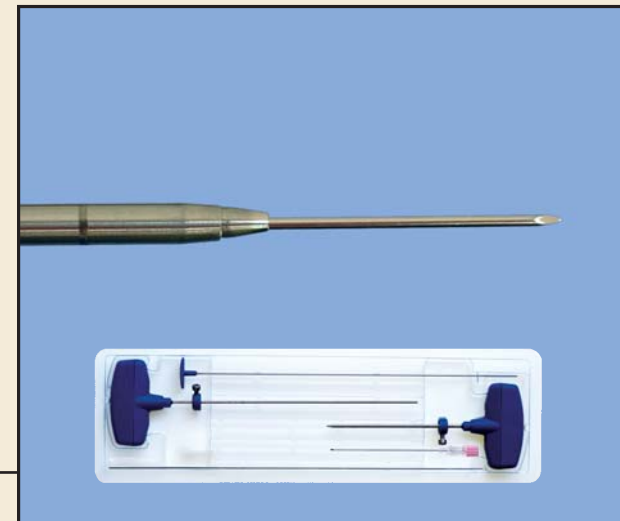

Coaxial Bone Biopsy Kit is mostly used to biopsy deep bones - Its design preserves soft tissue

Hard Bone Biopsy Kit: KBP

- ◆ Anesthesia Needle
- ◆ Stiffening Stylet
- ◆ Laurane™ Perforating Cannula Drill
- ◆ Biopsy Needle With Sharp Trepined Tip
- ◆ Needle Obturator

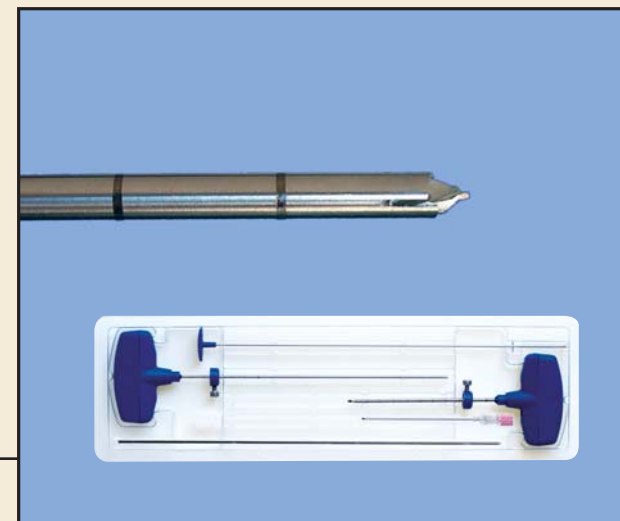

Hard Bone Biopsy Kit is used to biopsy hard bone - Its unique design* allows efficient penetration in hard bones.